

FUNCTION OF ROOT ELEVATORS IN SURGERY

To Loosen Tooth or Root From Bony Socket Prior to Placement of The Extraction Forceps.

HANDLES FOR ROOT ELEVATORS POINTS A, B, C, D, E, F, G, H

Thumb Handle



Shortblast Handle



Half Gold Handle



Round With Edge Handle



Thom with Flat Handle



Eye Style Handle



Eye Style Handle



Grip Style Handle



ROOT ELEVATORS

FUNCTION OF ROOT ELEVATORS IN SURGERY

To Loosen Tooth or Root From Bony Socket Prior to Placement of The Extraction Forceps.

Dental Root Elevators

To look at one of these instruments you might mistake it for a narrow screwdriver.

That's because just like one an elevator has a handle and then a specially designed "blade" or tip portion.

How they're used.

The tip of an elevator is designed to be wedged into the ligament space between a tooth and its surrounding bone (as shown in our animation).

As the elevator is forced into and twisted around in this space, the tooth is rocked around and pressed against the walls of its socket.

This helps to both expand the shape of the socket and separate the tooth from its ligament.

As this work is continued, the tooth gradually becomes more and more mobile.

In some cases, the dentist may be able to completely remove the tooth just using their elevator.

If not, they will choose a point where they feel they have accomplished as much as they can with it and will then switch to using extraction forceps.



■ Элеватор корневой экстракционный ■ Слеватор коренвий екстракційний
 ■ Elévateurs à racines ■ Elevadores para raíces
 ■ Wurzelheber ■ المصاعد الجذر

Article# IC5-001 LUXATOR ROOT ELEVATORS SET OF 6 PCS

FUNCTION: To Loosen Tooth or Root From Bony Socket Prior to Placement of The Extraction Forceps.



Article# IC5-002
Luxator
3mm
straight



Article# IC5-003
Luxator
4mm
straight



Article# IC5-004
Luxator
5mm
straight



Article# IC5-005
Luxator
3mm
curved



Article# IC5-006
Luxator
4mm
curved



Article# IC5-007
Luxator
5mm
curved



Article# IC5-008 LUXATOR ROOT ELEVATORS SET OF 8 PCS

FUNCTION: To Loosen Tooth or Root From Bony Socket Prior to Placement of The Extraction Forceps.



Article# IC5-009
Luxator
2mm
straight



Article# IC5-010
Luxator
3mm
straight



Article# IC5-011
Luxator
4mm
straight



Article# IC5-012
Luxator
5mm
straight



Article# IC5-013
Luxator
2mm
curved



Article# IC5-014
Luxator
3mm
curved



Article# IC5-015
Luxator
4mm
curved



Article# IC5-016
Luxator
5mm
curved



Article# IC5-017 LUXATOR ROOT ELEVATORS SET OF 6 PCS

FUNCTION: To Loosen Tooth or Root From Bony Socket Prior to Placement of The Extraction Forceps.



Article# IC5-018
Luxator
2mm
straight



Article# IC5-019
Luxator
3mm
straight



Article# IC5-020
Luxator
4mm
straight



Article# IC5-021
Luxator
2mm
curved



Article# IC5-022
Luxator
3mm
curved



Article# IC5-023
Luxator
4mm
curved



- Элеватор корневой экстракционный
- Elévateurs à racines
- Wurzelheber
- Слеватор корневий екстракційний
- Elevadores para raíces
- المساعد الجذر

Function of Root Elevators:

To loosen tooth or root from bony socket prior to placement of the extraction forceps.

ROOT ELEVATORS



Article# IC5-024
Luxating Elevators
1.5mm, Straight



Article# IC5-025
Luxating Elevators
2mm, Straight



Article# IC5-026
Luxating Elevators
2mm, Curved



Article# IC5-027
Luxating Elevators
3mm, Straight



Article# IC5-028
Luxating Elevators
3mm, Curved



Article# IC5-029
Luxating Elevators
3mm, Double Curved



Article# IC5-030
Luxating Elevators
3mm, Curved Up



Article# IC5-031
Luxating Elevators
4mm, Straight



Article# IC5-032
Luxating Elevators
4mm, Curved




Article# IC5-033
Luxating Elevators
5mm, Straight



Article# IC5-034
Luxating Elevators
5mm, Curved


This section includes the Luxating Elevator set with Stainless Steel Handle. Long & light weight handle design provide excellent smooth grip with controlled force transfer. The working procedure is similar as mentioned in the previous section.

- | | |
|------------------------------------|-------------------------------------|
| ■ Элеватор корневой экстракционный | ■ Elevator corneviu ekstraktsionnyi |
| ■ Elévateurs à racines | ■ Elevadores para raíces |
| ■ Wurzelheber | ■ المصاعد الجذر |

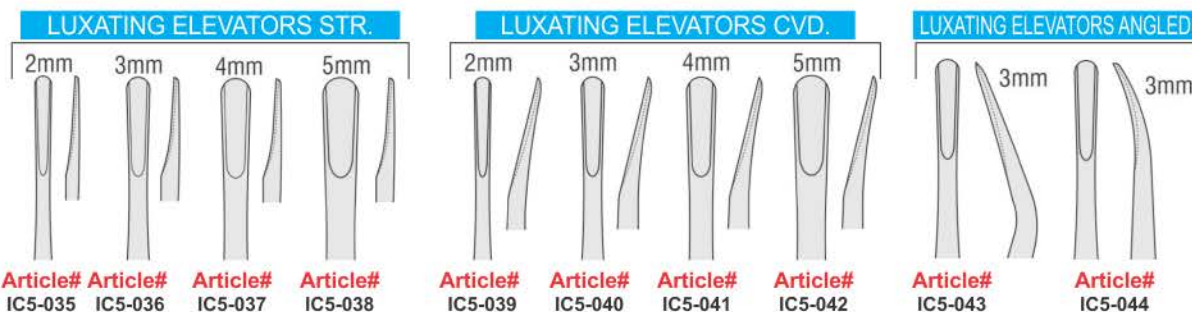

Function of Root Elevators:

To loosen tooth or root from bony socket prior to placement of the extraction forceps.

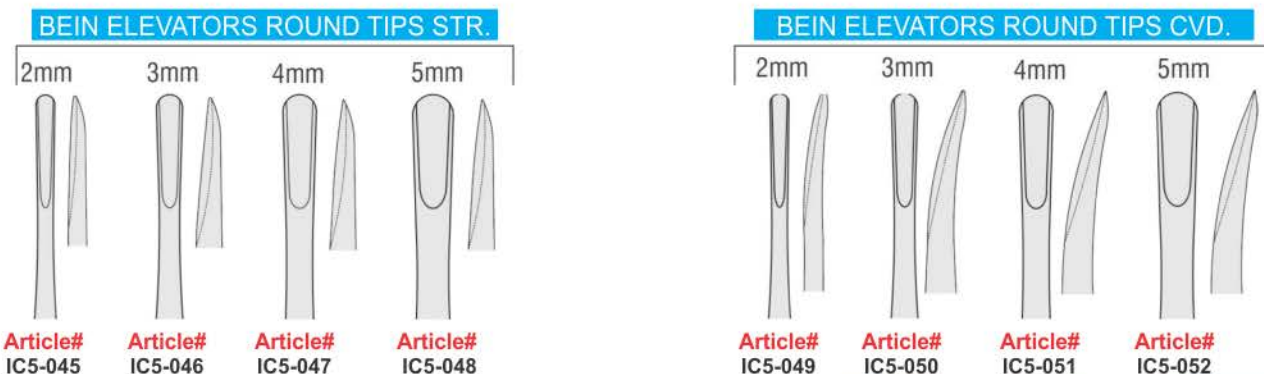

Upper root elevator with helicoidal tip

lower roots elevator


universal mini-elevator for fragments

Function of Root Elevators:

To loosen tooth or root from bony socket prior to placement of the extraction forceps.

■ Элеватор корневой экстракционный ■ Elevator corneviu extracțiunii
 ■ Elévateurs à racines ■ Elevadores para raíces
 ■ Wurzelheber ■ المصاعد الجذر

BEIN ELEVATORS POINTED TIPS STR.

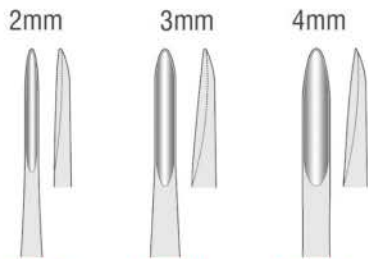


Article#
IC5-053

Article#
IC5-054

Article#
IC5-055



BEIN ELEVATORS POINTED TIPS CVD.

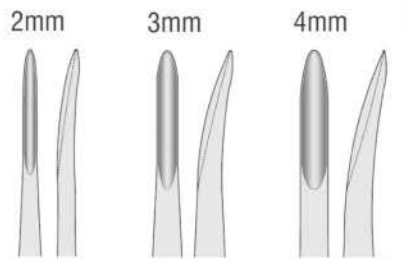

Article#
IC5-056

Article#
IC5-057

Article#
IC5-058

Function of Root Elevators:

To loosen tooth or root from bony socket prior to placement of the extraction forceps.

Straight



Curved



ROOT ELEVATORS



Article#
IC5-059
Bein Root
Elevator

Article#
IC5-060
Bein Root
Elevator

Article#
IC5-061
Bein Root
Elevator

Article#
IC5-062
Bein Root
Elevator




Article#
IC5-063
Bein Root
Elevator

Article#
IC5-064
Bein Root
Elevator

Article#
IC5-065
Bein Root
Elevator

Article#
IC5-066
Bein Root
Elevator

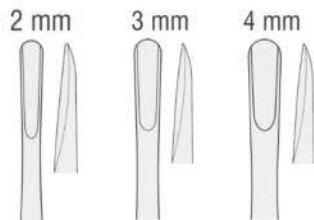


Article#
IC5-068
Bein
3mm

Article#
IC5-069
Bein
4mm



Article#
IC5-067
Bein 2mm

Bein Root Elevator

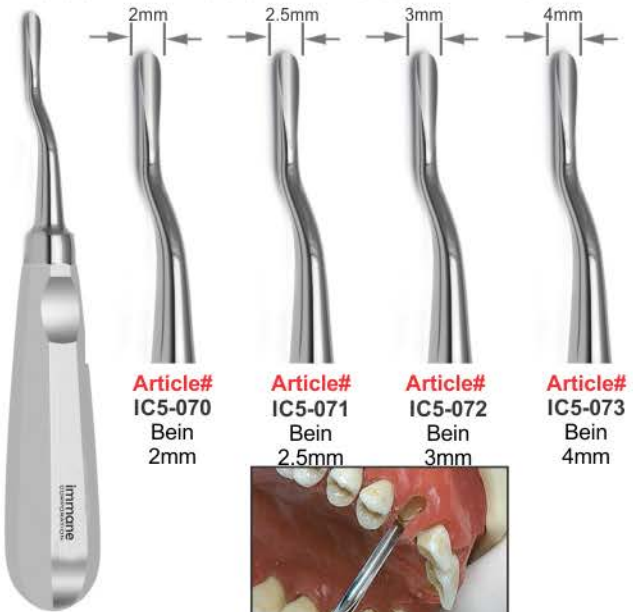

Article#
IC5-070
Bein
2mm

Article#
IC5-071
Bein
2.5mm

Article#
IC5-072
Bein
3mm

Article#
IC5-073
Bein
4mm


Bein Root Elevator



Элеватор корневой экстракционный Элеватор корневий екстракційний
Elévateurs à racines Elevadores para raíces


Function of Root Elevators:

To loosen tooth or root from bony socket prior to placement of the extraction forceps.

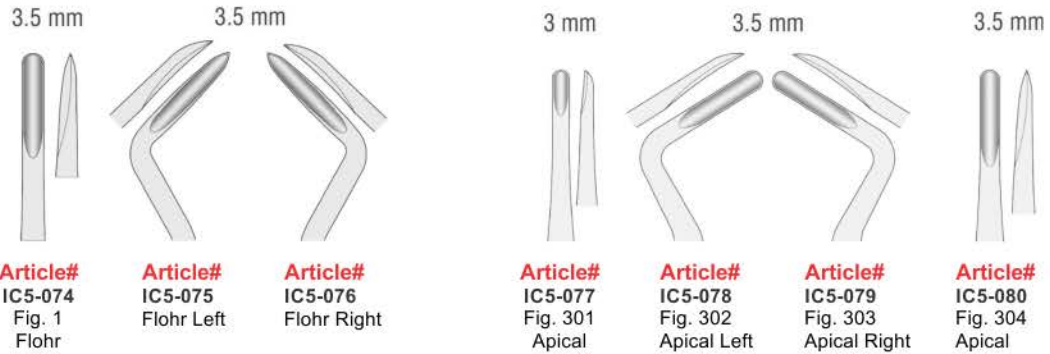

Function of Root Elevators:

To loosen tooth or root from bony socket prior to placement of the extraction forceps.

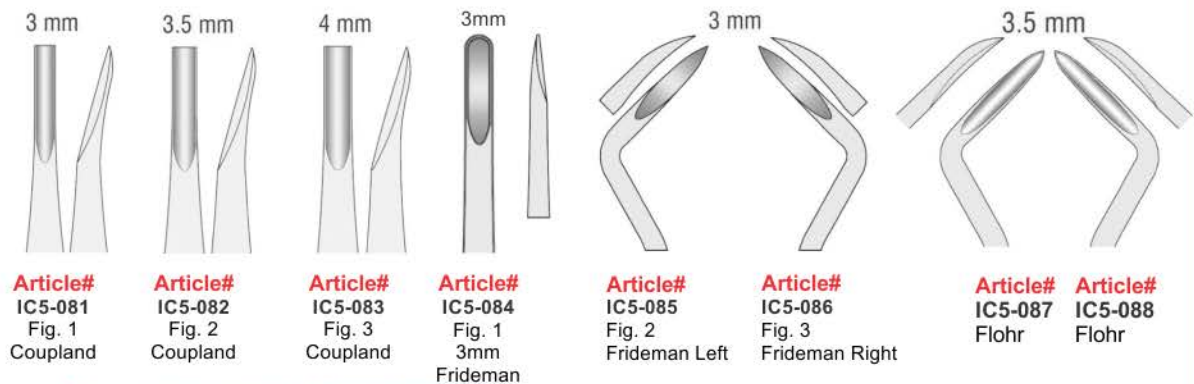

Apical Elevators

Apical root elevators, right, left
For upper and lower root tips
Width of tip 3 mm
Smooth Surface Handle
Surgical Stainless Steel
Sterilized in all Commercial Sterilizers



- | | |
|----------------------------------|---------------------------------|
| Элеватор корневой экстракционный | Слеватор корневий екстракційний |
| Elévateurs à racines | Elevadores para raíces |
| Wurzelheber | المصاعد الجذر |

Function of Root Elevators:

To loosen tooth or root from bony socket prior to placement of the extraction forceps.



Article# IC5-092
5
Schmeckebier
Apexo

Article# IC5-093
71 Miller
Apexo

Article# IC5-094
72 Miller
Apexo

Article# IC5-095
6 Potts
For impacted upper 3rd Molars
without bony obstruction.

Article# IC5-096
7 Potts

Article# IC5-097
73 Miller
Apexo

Article# IC5-098
74 Miller
Apexo

Article# IC5-099
19
Shorter leverage for lower 1st and 2nd
Molar roots. Also useful for premolar
and molar deciduous teeth.

Article# IC5-100
20



Article# IC5-091
4 Schmeckebier Apexo
Schmeckebier and Miller
designs used for displacement
of broken or deep-seated root tips.



Article# IC5-101
38
Serrated



Article# IC5-102
39
Serrated



Article# IC5-103
190



Article# IC5-104
191



Article# IC5-105
15
Woodward
Lower Molar roots



Article# IC5-106
16
Woodward
Lower Molar roots

For removal of root fragments

Элеватор корневой экстракционный Слеватор корневий екстракційний
Elévateurs à racines Elevadores para raíces
Wurzelheber المصاعد الجذر



Cryer root or fragment elevator



Seldin root or fragment elevator

Function of Root Elevators:

To loosen tooth or root from bony socket prior to placement of the extraction forceps.

Medium Cryer Root Elevators set
Left & Right, with medium working points
Working Tips are smooth
Smooth Surface Handle
Surgical Stainless Steel
Sterilized in all Commercial Sterilizers



Seldin root or fragment elevator



Cryer root or fragment elevator

Large Cryer Root Elevators set
Left & Right, with large working points
Working Tips are smooth
Smooth Surface Handle
Surgical Stainless Steel
Sterilized in all Commercial Sterilizers



Article# IC5-113
CRYER
44 Small

- | | |
|----------------------------------|---------------------------------|
| Элеватор корневой экстракционный | Слеватор корневий екстракційний |
| Elévateurs à racines | Elevadores para raíces |
| Wurzelheber | المساعد الجذر |

Function of Root Elevators:

To loosen tooth or root from bony socket prior to placement of the extraction forceps.

LENDO-LEVIAN ELEVATORS STR.



Article#
IC5-123

Article#
IC5-124

Article#
IC5-125

LENDO-LEVIAN ELEVATORS CVD.



Article#
IC5-126
Fig. 1

Article#
IC5-127
Fig. 2

Article#
IC5-128
Fig. 3

These root elevators are used for releasing the ligaments around the neck of the tooth, and to facilitate the extraction. The working end serration can be very helpful during the initial phase of extraction.

Elevator according to Lindo-Levien
Working part straight with serration
Width 2.2mm 3.0mm or 4.0mm
Smooth Surface Handle
Surgical Stainless Steel
Sterilized in all Commercial Sterilizers



Lindo-Levien Elevators

Upper root elevator with helicoidal tip

lower roots elevator



Article# IC5-129
Lindo-Levien
Elevators
3mm, Straight

Article# IC5-130
Lindo-Levien
Elevators
3mm, Curved

Article# IC5-131
Lindo-Levien
Elevators
4mm, Straight

Article# IC5-132
Lindo-Levien
Elevators
4mm, Curved



Article#
IC5-133
Heidbrink
Left

Article#
IC5-134
Heidbrink
Right

HEIDBRINK
Root Tip Pick

Article# IC5-135
Heidbrink Straight

- | | |
|------------------------------------|-------------------------------------|
| ■ Элеватор корневой экстракционный | ■ Elevator cornevei ekstraktsionnyi |
| ■ Elévateurs à racines | ■ Elevadores para raíces |
| ■ Wurzelheber | ■ المصاعد الجذر |

Coupland Elevators

Function of Root Elevators:

To loosen tooth or root from bony socket prior to placement of the extraction forceps.

Coupland root elevators set 3mm, 3.5mm & 4mm tip width
Sand Blasted no slip handle Surgical Stainless Steel
Sterilized in all Commercial Sterilizers



Article# IC5-136



Article# IC5-137



Article# IC5-138



Article# IC5-139



Article# IC5-140



Article# IC5-141



Article# IC5-142

Seldin Elevators

Seldin root elevators Smooth Surface Handle
Surgical Stainless Steel Sterilized in
all Commercial Sterilizers

B-Cogswell & Apexo Elevators

Smooth Surface Handle Surgical
Stainless Steel Sterilized in
all Commercial Sterilizers

Potts Elevators
These root elevators are used for releasing the
ligaments around the neck of the tooth, and to
facilitate the extraction. The working end are
smooth These are used for impacted upper
3rd Molars without bony obstruction.

Elevator according to Potts
Working Tips are smooth
Smooth Surface Handle
Surgical Stainless Steel
Sterilized in all Commercial Sterilizers

Wisdom Roots Elevator

34 Str. Seldin Elevator, Upper and
lower anterior roots



Article#
IC5-143

3S Seldin Elevator, Upper & lower
anterior roots, sharp cone edge



Article#
IC5-144

2S Seldin Elevator, Lower 3rd Molars



Article#
IC5-145

301 Apexo Elevator, for displacement of deeply seated roots



Article#
IC5-146

B Cogswell Elevator



Article#
IC5-147

Article# IC5-148
Wisdom for 28
and 48




Article# IC5-149
Wisdom for 18
and 38



■ Элеватор корневой экстракционный ■ Слеватор кореневий екстракційний
 ■ Elévateurs à racines ■ Elevadores para raíces
 ■ Wurzelheber ■ المصاعد الجذر

ROOT ELEVATORS

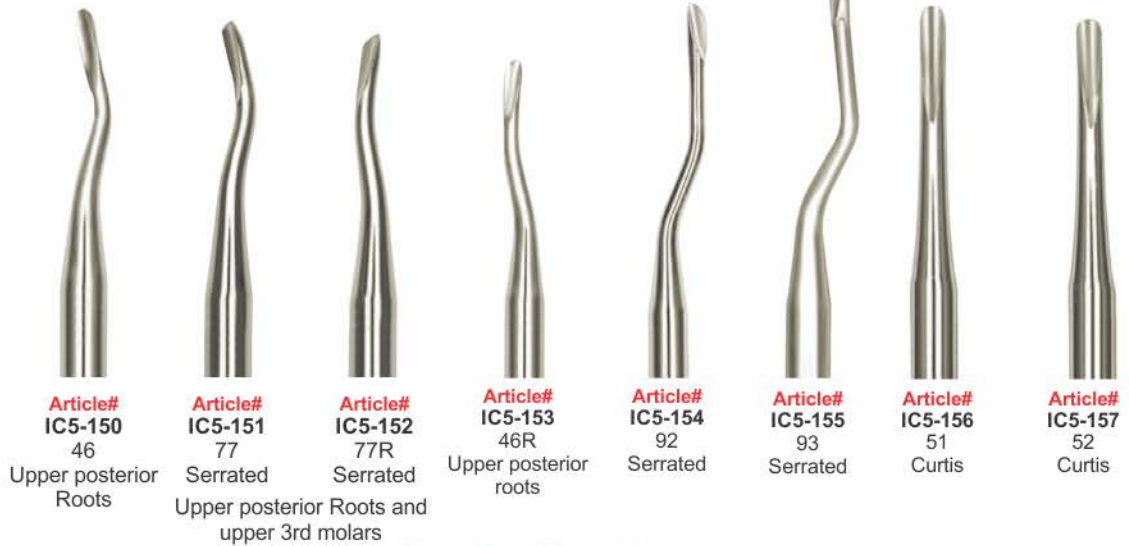

Article# IC5-150 46 Upper posterior Roots
Article# IC5-151 77 Serrated Upper posterior Roots and upper 3rd molars
Article# IC5-152 77R Serrated
Article# IC5-153 46R Upper posterior roots
Article# IC5-154 92 Serrated
Article# IC5-155 93 Serrated
Article# IC5-156 51 Curtis
Article# IC5-157 52 Curtis

Upper and lower molar roots **Function of Root Elevators:**
To loosen tooth or root from bony socket prior to placement of the extraction forceps.



Article# IC5-159 34 Upper and lower anterior roots
Article# IC5-160 2S Lower 3rd Molars
Article# IC5-161 3S Upper and lower anterior roots
Article# IC5-162 8 Crane Pick
Article# IC5-163 44C Bradley For displacement of roots and 3rd Molars
Article# IC5-164 A Cogswell
Article# IC5-165 B Cogswell

Article# IC5-158 34S
 Thinner tip and more narrow than the 34 shank

Function of Root Elevators:
To loosen tooth or root from bony socket prior to placement of the extraction forceps.

Seldin root or fragment elevator
 3.5mm



Cryer root or fragment elevator



Article# IC5-166 Seldin
Article# IC5-167 11 Stout 3rd Molars
Article# IC5-168 11A Stout 3rd Molars

- | | |
|----------------------------------|---------------------------------|
| Элеватор корневой экстракционный | Слеватор корневий екстракційний |
| Elévateurs à racines | Elevadores para raíces |
| Wurzelheber | المساعد الجذر |

Function of Root Elevators:

To loosen tooth or root from bony socket prior to placement of the extraction forceps.

Laxutor Root Elevators Set of 9 Pcs



Article# IC5-169
2mm Str.



This section includes a comprehensive set of Posterior Flared Tip Elevators. These elevators used to loosen the tooth from the periodontal ligament and ease extraction, particularly used for Upper Posterior Roots and Upper 3rd Molar.

E46 & E46R: Mesial & Distal with 2.75mm, Smooth working tips

E77 & E77R: Mesial & Distal with 3.75mm, Serrated working tips

Surgical Stainless Steel Hollow handle, very light weight Completely Sterilizable



Posterior Flared Tip Elevators



Cryer root or fragment elevator



Article# IC5-175
3mm Angle



Article# IC5-176
3mm Str.



Article# IC5-177
5mm Str.



Luxating Elevator With Periotome Tip



Article# IC5-178
46R Elavator, Smooth,
Mesial, 2.75mm



Article# IC5-179
46 Elavator, Smooth,
Distal, 2.75mm



Article# IC5-180
77R Elavator, Serrated,
Mesial, 3.75mm



Article# IC5-181
77 Elavator, Serrated,
Distal, 3.75mm



Article# IC5-182
3mm Luxating
Elevator / Periotome



Article# IC5-183
4mm
Luxating Elevator
Periotome

Article# IC5-184
5mm Luxating
Elevator
Periotome



- | | |
|----------------------------------|---------------------------------|
| Элеватор корневой экстракционный | Слеватор корневий екстракційний |
| Elévateurs à racines | Elevadores para raíces |
| Wurzelheber | المساعد الجذر |

Function of Root Elevators:

To loosen tooth or root from bony socket prior to placement of the extraction forceps.
Kopp Roots Elevator



Article# IC5-185
Fig. 1
Kopp 3.6mm Upper Roots



Article# IC5-186
Fig. 2
Kopp 3.6mm Upper Roots



Article# IC5-187
Fig. 3
Kopp 4.2mm Upper Roots



Article# IC5-188
Fig. 4
Kopp 4.2mm Upper Roots



Article# IC5-189
Fig. 1 Kopp 3.6mm Lower Roots



Article# IC5-190
Fig. 2 Kopp 3.6mm Lower Roots



Article# IC5-191
Fig. 3 Kopp 4.2mm Lower Roots



Article# IC5-192
Fig. 4 Kopp 4.2mm Lower Roots



lower roots elevator



Cryer root or fragment elevator

Apexo Root Elevators



Article# IC5-201
Bein



Article# IC5-193
Bein



Article# IC5-194
Bein



Article# IC5-195
Bein



Article# IC5-202
Gartner



Article# IC5-203
Gartner



Article# IC5-204
Gartner



Article# IC5-196
Fig. 301
Apexo 2.5mm



Article# IC5-197
Fig. 302
Apexo 2.5mm



Article# IC5-198
Fig. 303
Apexo 2.5mm



Article# IC5-199
Fig. 301W
Seldin 3mm



Article# IC5-200
Fig. 304W
Seldin 3mm



universal mini-elevator for fragments

- | | |
|------------------------------------|----------------------------------|
| ■ Элеватор корневой экстракционный | ■ Elevator corneviu extracționii |
| ■ Elévateurs à racines | ■ Elevadores para raíces |
| ■ Wurzelheber | ■ المصاعد الجذر |

Function of Root Elevators:

To loosen tooth or root from bony socket prior to placement of the extraction forceps.

Medan Roots Elevator

Hylin Roots Elevator



- | | | | | | |
|----------------------------|----------------------------|----------------------------|----------------------------|----------------------------|----------------------------|
| Article#
IC5-205 | Article#
IC5-206 | Article#
IC5-207 | Article#
IC5-208 | Article#
IC5-209 | Article#
IC5-210 |
| ein 2.5mm | Bein 3mm
Flat | Bein 3mm
Round | Bein 3.5mm | Flohr 3.5mm
(Left) | Flohr 3.5mm
(Right) |



- | | | |
|----------------------------|----------------------------|----------------------------|
| Article#
IC5-211 | Article#
IC5-212 | Article#
IC5-213 |
| Hylin
4mm
Str. | Hylin
3.5mm
Left | Hylin
3.5mm
Right |



Cryer root or fragment elevator

Medan Roots Elevator

- Article#** IC5-214
Hylin
3mm
Str.

Overall length
154mm

HEIDBRINK



- | | | | | | | | | | | |
|----------------------------|----------------------------|----------------------------|----------------------------|----------------------------|---|----------------------------|----------------------------|----------------------------|----------------------------|----------------------------|
| Article#
IC5-216 | Article#
IC5-217 | Article#
IC5-218 | Article#
IC5-219 | Article#
IC5-220 | Article#
IC5-221 | Article#
IC5-222 | Article#
IC5-223 | Article#
IC5-224 | Article#
IC5-225 | Article#
IC5-226 |
| 303 | 301W | 304W | 1 | 2 | 80 | 81 | 81A | 81E | 40 | 41 |
| Apexo | Seldin | Seldin | Heidbrink
Straight | Heidbrink
Straight | Upper and lower anterior teeth
and roots | | | Heidbrink
Straight | | |

Upper and lower anterior teeth and roots
Upper anterior: For displacement of anterior roots
and small teeth



- Article#** IC5-215
301 Apexo

Upper root elevator with helicoidal tip
For displacement of deeply seated roots

lower roots elevator

universal mini-elevator for fragments

- | | |
|----------------------------------|---------------------------------|
| Элеватор корневой экстракционный | Слеватор корневий екстракційний |
| Elévateurs à racines | Elevadores para raíces |
| Wurzelheber | المساعد الجذر |

Function of Root Elevators:

To loosen tooth or root from bony socket prior to placement of the extraction forceps.



Root Elevators



FUNCTION: To Loosen Tooth or Root From Bony Socket Prior to Placement of The Extraction Forceps.

Mesial & Distal Root Elevators



Root Elevators



- | | |
|----------------------------------|---------------------------------|
| Элеватор корневой экстракционный | Слеватор корневий екстракційний |
| Elévateurs à racines | Elevadores para raíces |
| Wurzelheber | المساعد الجذر |

Function of Root Elevators:

To loosen tooth or root from bony socket prior to placement of the extraction forceps.

Bernard Root Elevators are particularly suitable for the excision of fractured roots or root fragments. They represent a combination of root tip elevators and luxating elevators. The inner side is concave and the lateral edges are for cutting. This allow for rotatory movements, and when applying slight pressure, the ligaments can be separated and the remaining root fragments can be loosened.

This section includes a set of five Bernard Root Elevators, including Left, Right, Straight, Mesial & Distal form. All elevators are made from, Surgical Stainless Steel Hallow handle, light weight Completely Sterilizable

Bernard Root elevators also know as Syndesmotomes. They can also be use as Periosteum or to elevate Small Roots



Bernard Elevators -Syndesmotomes

Article# IC5-255 BAYONET ELEVATORS SET OF 7 PCS



- | | |
|----------------------------------|---------------------------------|
| Элеватор корневой экстракционный | Слеватор корневий екстракційний |
| Elévateurs à racines | Elevadores para raices |
| Wurzelheber | المساعد الجذر |

Function of Root Elevators:

To loosen tooth or root from bony socket prior to placement of the extraction forceps.

Article# IC5-266

LUXATORS SET OF 9 PCS



Article# IC5-276

LUXATORS SET OF 7 PCS



Элеватор корневой экстракционный Слеватор корневий екстракційний
Elévateurs à racines Elevadores para raíces
Wurzelheber المصاعد الجذر

ST22284

LUXATORS SET OF 4 PCS


Directa

Description:


The Luxator P-series offers good tactility and reach with the penholder grip. These instruments are thinner and sharper than the original Luxator instruments. Initiating the extraction process with the P-series range will simplify your work. Luxator P-series kit Contains: 4 different blade shapes: ST22415 Luxator P1, straight blade, dark green - Straight universal periosteum, ideal for anteriors and premolars. ST22416 Luxator P2, curved blade, light blue - A slightly curved periosteum for improved ergonomics. ST22417 Luxator P3, contra angle blade, dark blue - Distal use, specially designed for use on posteriors. ST22418 Luxator P4, dual edge blade, lilac - For fractured root tips and very decayed teeth




Article#
IC5-285

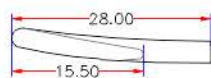

Article#
IC5-286

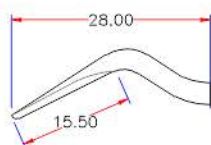


Article#
IC5-287

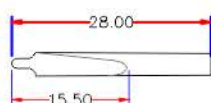


Article#
IC5-288


Article# IC5-289 Luxating Elevators Set of 5 pcs

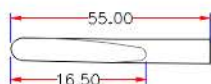


Article# IC5-290 Luxating Elevator #2, 3mm Straight

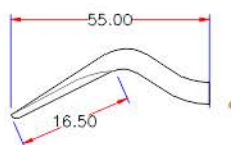


Article# IC5-291 Luxating Elevator #7, 3mm Double Curved

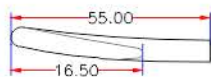

Article# IC5-292 Luxating Elevator #4, 3mm Curved

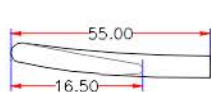

Article# IC5-293 Luxating Elevator #5, 5mm Curved

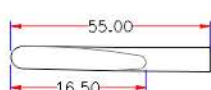





Article# IC5-294 Luxating Elevator #3, 5mm Straight

- | | |
|------------------------------------|-----------------------------------|
| ■ Элеватор корневой экстракционный | ■ Слеватор корневий екстракційний |
| ■ Elévateurs à racines | ■ Elevadores para raíces |
| ■ Wurzelheber | ■ المصاعد الجذر |

Function of Root Elevators:

To loosen tooth or root from bony socket prior to placement of the extraction forceps.

LUXATORS SET OF 10 PCS MICRO SERRATED

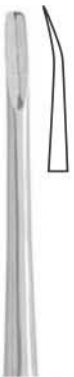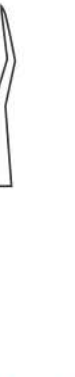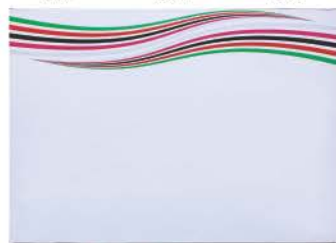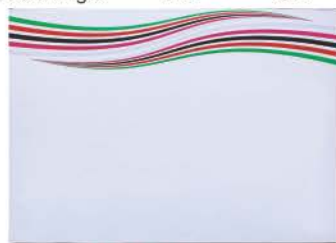											
	Article# IC5-296 Laxutor 1.5mm Str.	Article# IC5-297 Laxutor 2mm Str.	Article# IC5-298 Laxutor 3mm Str.	Article# IC5-299 Laxutor 5mm Str.	Article# IC5-300 Laxutor 3mm Cvd.	Article# IC5-301 Laxutor 5mm Cvd.	Article# IC5-302 Laxutor 3mm Upper Angle	Article# IC5-303 Laxutor 3mm Double Angle	Article# IC5-304 Laxutor 3mm Str.	Article# IC5-305 Laxutor 5mm Str.	
											
	Article# IC5-295 LUXATORS SET OF 10 PCS MICRO SERRATED				Article# IC5-306 BEIN ROOT ELEVATORS SET OF 6 PCS			Article# IC5-306 BEIN ROOT ELEVATORS SET OF 6 PCS			

					
Article# IC5-307 Bein 3S 3mm	Article# IC5-308 Bein 3IC 2mm	Article# IC5-309 Bein 3C 3mm	Article# IC5-310 Bein 4C 4mm	Article# IC5-311 Bein 5C 3mm	Article# IC5-312 Bein 5CA 3mm

- | | |
|------------------------------------|-------------------------------------|
| ■ Элеватор корневой экстракционный | ■ Elevator corneviu ekstraktsionnyi |
| ■ Elévateurs à racines | ■ Elevadores para raíces |
| ■ Wurzelheber | ■ المصاعد الجذر |

Function of Root Elevators:

To loosen tooth or root from bony socket prior to placement of the extraction forceps.



Article# IC5-313
301
Apexo



Article# IC5-314
301
Apexo
Ergo Dynamic



Article# IC5-315
#77R Elevator
Serrated
Ergo Dynamic



Article# IC5-316
#B Cogswell
Elevator
Ergo Dynamic



Article# IC5-317
#34S
Seldin Elevator
Straight
Ergo Dynamic



Article# IC5-318
Luxating Elevator
Straight
3mm
Ergo Dynamic



Article# IC5-319
Luxating Elevator
Curved
3mm
Ergo Dynamic

ERGO DYNAMIC ELEVATORS OFFER AN ALTERNATIVE TO STANDARD SIZE ELEVATOR HANDLES.



Article# IC5-320
Luxating Elevator
Straight
5mm
Ergo Dynamic



Article# IC5-321
Luxating Elevator
Curved
5mm
Ergo Dynamic

LUXATOR ROOT ELEVATORS SET OF 6 PCS



Article# IC5-322



Article# IC5-323



Article# IC5-324



Article# IC5-325



Article# IC5-326



Article# IC5-327



Article# IC5-328
Root Elevator Set of 6 Pcs

- | | |
|----------------------------------|---------------------------------|
| Элеватор корневой экстракционный | Слеватор корневий екстракційний |
| Elévateurs à racines | Elevadores para raíces |
| Wurzelheber | المساعد الجذر |

Function of Root Elevators:

To loosen tooth or root from bony socket prior to placement of the extraction forceps.

FUNCTION: To Loosen Tooth or Root From Bony Socket Prior to Placement of The Extraction Forceps.

W. Barry Roots Elevator



W. Barry
Fig. 1
(125mm)

Article# IC5-329
W. Barry
Fig. 1

Article# IC5-330
W. Barry
Fig. 2

Article# IC5-331
W. Barry
Fig. 3

Article# IC5-332
W. Barry
Fig. 4



Winter Roots Elevator



Winter
Fig. 1L
(125mm)



Article#
IC5-333
Winter
Fig. 1L



Article#
IC5-334
Winter
Fig. 1R



Article#
IC5-335
Winter
Fig. 11L



Article#
IC5-336
Winter
Fig. 11R



Article#
IC5-337
Winter
Fig. 12L



Article#
IC5-338
Winter
Fig. 12R

Pott Roots Elevator



Article# IC5-339
Winter
Fig. 1L



Article# IC5-340
Winter
Fig. 1R



Article# IC5-341
Lecluse
bayonet



Article# IC5-342
Pott
Left



Article# IC5-343
Pott
Right

- | | |
|------------------------------------|-----------------------------------|
| ■ Элеватор корневой экстракционный | ■ Слеватор корневий екстракційний |
| ■ Elévateurs à racines | ■ Elevadores para raíces |
| ■ Wurzelheber | ■ المصاعد الجذر |

Function of Root Elevators:

To loosen tooth or root from bony socket prior to placement of the extraction forceps.

FUNCTION: To Loosen Tooth or Root From Bony Socket Prior to Placement of The Extraction Forceps.

Heidbrink Roots Elevator

Syndesmotomes Roots Elevator



Periotome Roots Elevator

Syndesmotomes



- | | |
|----------------------------------|---------------------------------|
| Элеватор корневой экстракционный | Слеватор корневий екстракційний |
| Elévateurs à racines | Elevadores para raíces |
| Wurzelheber | المساعد الجذر |

Function of Root Elevators:

To loosen tooth or root from bony socket prior to placement of the extraction forceps.

HOWARD Slightly heavier design than Heidbrink picks.

HEIDBRINK



Article#
IC5-362
Heidbrink
Root Apex
Elevator
(Left)

Article#
IC5-361
Heidbrink
Root Apex
Elevator



Article# IC5-363
Heidbrink
Root Apex
Elevator (Right)



Article#
IC5-364
2
Howard



Article#
IC5-366
3
Howard



Article#
IC5-365
1
Howard



Article#
IC5-367
2
Heidbrink



Article#
IC5-369
3
Heidbrink

Article#
IC5-368
1
Heidbrink

FUNCTION: To loosen small root fragments from bony socket

Heidbrink Roots Elevator

Warwick-James



Article#
IC5-370/1
Heidbrink
Fig.1

Article#
IC5-370/2
Heidbrink
Fig.2

Article#
IC5-370/3
Heidbrink
Fig.3

Heidbrink
Fig.1



Article#
IC5-372
Fig. 2

Article#
IC5-371
Fig. 1

Article#
IC5-373
Fig. 3



Article#
IC5-374
Fig. 1
Chompret



Article#
IC5-375
Fig. 2
Chompret



Article#
IC5-376
Fig. 3
Chompret



Article#
IC5-377
Fig. 4
Chompret



Chompret syndesmotome, sickle pattern

- | | |
|------------------------------------|-----------------------------------|
| ■ Элеватор корневой экстракционный | ■ Слеватор корневий екстракційний |
| ■ Elévateurs à racines | ■ Elevadores para raíces |
| ■ Wurzelheber | ■ المصاعد الجذر |

Function of Root Elevators:

To loosen tooth or root from bony socket prior to placement of the extraction forceps.



Article#
IC5-378
Heidbrink
Left



Article# **IC5-380**
Heidbrink
Straight



Article#
IC5-379
Heidbrink
Right



Article#
IC5-381
Warwick-
James
Straight



Article#
IC5-382
Warwick-
James
Left



Article#
IC5-383
Warwick-
James
Right

FUNCTION: To Loosen Tooth or Root From Bony Socket Prior to Placement of The Extraction Forceps.



Coupland



Article#
IC5-384
Fig. 1
190/15



Article#
IC5-385
Fig. 2
190/16



Article#
IC5-386
Fig. 3
190/17



universal mini-elevator for fragments



Article#
IC5-387
192/1



Article#
IC5-388
192/2



Article#
IC5-389
192/3



Article#
IC5-390
Fig. 1



Article#
IC5-391
Fig. 2



Article#
IC5-392
Fig. 3

- | | |
|----------------------------------|---------------------------------|
| Элеватор корневой экстракционный | Слеватор корневий екстракційний |
| Elévateurs à racines | Elevadores para raíces |
| Wurzelheber | المساعد الجذر |

Function of Root Elevators:

To loosen tooth or root from bony socket prior to placement of the extraction forceps.

PDL ELEVATORS



- | | |
|----------------------------------|---------------------------------|
| Элеватор корневой экстракционный | Слеватор корневий екстракційний |
| Elévateurs à racines | Elevadores para raíces |
| Wurzelheber | المساعد الجذر |

Function of Root Elevators:

To loosen tooth or root from bony socket prior to placement of the extraction forceps.




Article# IC5-424
Root Elevator Set of 7 Pcs

- | | | | | | | |
|--|--|--|--|--|---|---|
| Article#
IC5-417
Bein
DL 32
for upper roots
round point. | Article#
IC5-418
Bein
DL 33
for upper roots
round point. | Article#
IC5-419
Bein
DL 34
for upper roots
round point. | Article#
IC5-420
Flohr
Fig. 3
for lower roots | Article#
IC5-421
Flohr
Fig. 4
for lower roots | Article#
IC5-422
Cryer-White
for lower roots,
wide pattern | Article#
IC5-423
Cryer-White
for lower roots,
wide pattern |
|--|--|--|--|--|---|---|



- | | | | | | | |
|--|--|--|--|--|---|---|
| Article#
IC5-425
Bein
DL 32
for upper roots
round point. | Article#
IC5-426
Bein
DL 33
for upper roots
round point. | Article#
IC5-427
Bein
DL 34
for upper roots
round point. | Article#
IC5-428
Flohr
Fig. 3
for lower roots | Article#
IC5-429
Flohr
Fig. 4
for lower roots | Article#
IC5-430
Cryer-White
for lower roots,
wide pattern | Article#
IC5-431
Cryer-White
for lower roots,
wide pattern |
|--|--|--|--|--|---|---|

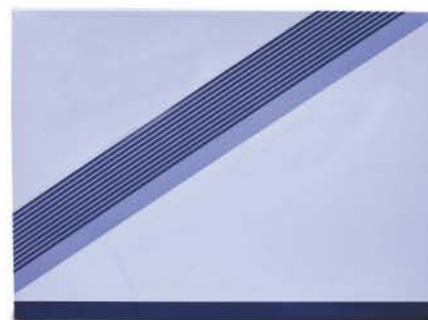

Article# IC5-432
Root Elevator Set of 7 Pcs

- | | |
|------------------------------------|-----------------------------------|
| ■ Элеватор корневой экстракционный | ■ Слеватор корневий екстракційний |
| ■ Elévateurs à racines | ■ Elevadores para raíces |
| ■ Wurzelheber | ■ المصاعد الجذر |

Function of Root Elevators:


To loosen tooth or root from bony socket prior to placement of the extraction forceps.



- | | | | | | | | | | |
|--|--|---|---|---|--|--|---|---|---|
| Article#
IC5-433
Root Elevator
DL 32
for upper roots
round point. | Article#
IC5-434
Root Elevator
DL 33
for upper roots
round point. | Article#
IC5-435
Bein
DL 34
for upper roots
round point. | Article#
IC5-436
Flohr
Fig. 3
for lower roots | Article#
IC5-437
Flohr
Fig. 4
for lower roots | Article#
IC5-438
Cryer-White
for lower roots,
wide pattern | Article#
IC5-439
Cryer-White
for lower roots,
wide pattern | Article#
IC5-440
Apical 301
3.5 mm | Article#
IC5-441
Apical 302
3.5 mm | Article#
IC5-442
Apical 303
3.5 mm |
|--|--|---|---|---|--|--|---|---|---|



Cryer root or fragment elevator

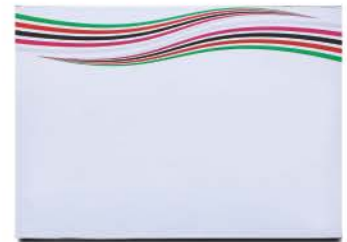

Article# IC5-443
Root Elevator Set of 7 Pcs



- | | | | | | | |
|---|---|---|---|--|---|--|
| Article#
IC5-444
Bein
DL 32
for upper roots
round point. | Article#
IC5-445
Bein
DL 33
for upper roots
round point. | Article#
IC5-446
Bein
DL 34
for upper roots
round point. | Article#
IC5-447
Cryer 27
(Left) | Article#
IC5-448
Cryer 28
(Right) | Article#
IC5-449
Flohr
3mm
(Left) | Article#
IC5-450
Flohr
3mm
(Right) |
|---|---|---|---|--|---|--|



Article# IC5-451
Root Elevator Set of 7 Pcs